



**MEDIC
THERAPEUTICS**
SPECIAL EDITION

VIBRATING FOOT MASSAGER




**MEDIC
THERAPEUTICS**
SPECIAL EDITION

User Manual

IMPORTANT

Read all instructions carefully before using this product. Retain this manual for future reference. Read the operating instructions carefully before use, especially the safety precautions.

SAFETY WARNINGS

- Always place product on a flat surface.
- Ensure that children do not use product without adult supervision.
- Do not use product if your hands are wet.
- Do not allow product to make contact with water.
- If the power cable is damaged, it should be replaced with a special cable by the manufacturer of authorized technical service.
- Unplug the device after use.
- As with any other physical exercise device, consult your doctor before you start using product.
- This device may be used by children aged eight and over, by people with reduced physical sensory or mental abilities, and by people lacking in experience and knowledge if they are appropriately supervised or given suitable instructions on the safe use of the device and understand the dangers its use entails. Children should not be allowed to play with the device. Children should not clean or do maintenance on device without supervision.
- This product is not intended for medical use. It is intended only for massages. The appliance should only be used with the power supply unit provided with the device.
- Consult your physician before using this product if:
 - You are pregnant
 - You have a pacemaker, an implant or other aids
 - You suffer from diabetes
 - You have a circulatory disorder
 - You have concerns about your health
 - The elderly and infirm should consult
- Because skin friction can be painful, it is recommended that the device not be directly applied to the skin. Stop using immediately if you feel discomfort or pain.
- Do not use this product at bedtime. The massager has stimulating effect and can delay sleep.

- To avoid a short circuit. Connect the device correctly and keep it away from metal (needles), waste and water.
- Use the device on a flat surface. If the device is not stable, it could tip over and cause injury.
- Remove the adapter from the power outlet when not in use.
- If you suspect a problem, turn off the device as soon as possible.
- When not in use for extended periods, store the device in its packaging in a safe, dry, cool place.
- Do not place any other objects on the massager.
- If the power cord is damaged, it must be replaced, but only by an authorized and qualified person.
- To avoid injury and irreparable damage, always sit while using the massager. To avoid dangerous injury, never stand on the unit or use it forcefully.
- Do not attempt to disassemble or repair the product. Repairs may only be carried out by suitably qualified technicians.


PRODUCT SPECIFICATIONS

Product Name:	VIBRATING FOOT MASSAGER
Input:	AC100-240V , 50/60Hz
Dimension:	40x36x13,2cm
Net weight:	3.9kg
Cross weight:	4.55kg

PRODUCT OVERVIEW




PRODUCT OVERVIEW

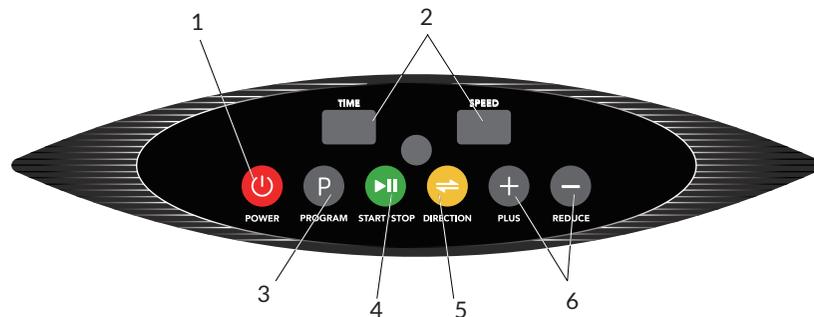

1. **Power:** Turn the Power On/Off
2. **LED Display Window:** displays the runtime and current speed
3. **Program:** three automatic programs
4. **Start/Stop:** start/stop the machine.
5. **Running Direction:** Adjust the machine running direction of the Platform while using it in manual mode
6. **+/-:** vibration speed adjustment.

NAMES AND FUNCTIONS OF THE CONTROL

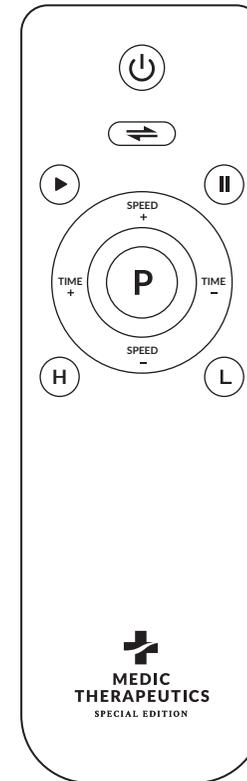

- POWER:**
Press to turn on/off
- RUNNING DIRECTION**
Adjust the machine running direction of the Platform while using it in manual mode
- STRAT:**
Press to strat the automatic setting
- STOP:**
Press to stop the automatic setting
- SPEED:**
Press to adjust the speed of the Platform while using it in manual mode
- TIME:**
Press to adjust the time of the duration before starting a new manula mode session
- H KEY:**
Starts and operates automatically at Level 12 speed. Note: the max speed of this foot massager is Level 15
- L KEY:**
Starts and operates automatically at Level 3 speed
- P KEY:**
The foot massager will automatically cycle through 3 preset programs:
 - P1- Low speed program (Levels 3-9)
 - P2- Middle speed program (Levels 6-12)
 - P3- High speed program (Levels 9-15)

Note: Each program operates for 10 minutes



Reception range of the remote control operation is 8 feet. Ensure that the red indicator light at the top of the remote control points towards the screen panel of the foot massager when using the remote control.

PRODUCT CARE

STORAGE

When not in use, place the product into the box and store in a dry, ventilated, cool area.

Do not twist the power cord.

Do not place product in direct sunlight or heat. Do not store where there is high humidity, in areas such as bathrooms.

CLEANING

Before cleaning, make sure to turn the power off and unplug the power cord.

Use a dry cloth to clean the surface of the product.

Do not use corrosive cleaning agents to clean the product.

TROUBLESHOOTING

If you have problems when using this product, refer to the following tips and possible solutions.

Status	Possible Cause	Solution
Product doesn't turn on	Not powered	Make sure power cord is plugged in
Product suddenly stops operating	<ol style="list-style-type: none">1. Power failure2. Timer exceeded 30 minutes3. Continuous use time is too long, overheat protection4. Max User Weight exerted on pedals has been exceeded. A Protection Action is initiated	<ol style="list-style-type: none">1. Check that the power cord is connected2. Press the Power button to turn the machine on3. Let the product cool down, then turn the power back on4. Do not stand on device and do not exceed the Max User Weight

LIMITED LIFETIME WARRANTY

Your Medic Therapeutics vibrating foot massager is backed by a limited lifetime manufacturer's warranty. Medic Therapeutics will repair or replace your device at any time should it fail due to a defect in material or workmanship, subject to the certain limitations.

This limited warranty does not cover any damage that results from unauthorized or improper use, service, or repair. Further, it does not cover damage caused by accident, impact, negligence, or normal wear and tear. Should you discover your Medic Therapeutics vibrating foot massager is not functioning properly, please send your device to our repair center for evaluation. If your product cannot be repaired or serviced, we will reserve the right to change it for a similar or newer model.

Please note that a flat rate of \$48.00 will be charged to cover a service evaluation fee and return shipping of your device. All warranty claims must be accompanied by a copy of your proof of purchase from an authorized retailer. Please send your device, proof of purchase, and a check or money order in the amount of \$48.00 made out to Medic Therapeutics to:

Address:

Medic Therapeutics Service Center
3069 Taft Street
Hollywood, FL 33021

Contact:

warranty@medictherapeutics.com