



FULL LEG THERMO AIR COMPRESSOR MASSAGER



User Manual

SAFETY PRECAUTIONS

If you have any of the following conditions, please consult your doctor prior to using:

1. Pregnant women or women after recently have given birth
2. Heart disease or implanted cardiac pacemaker
3. Fragile or damaged bones
4. Suffer from malignant tumors
5. Skin infections, irritation, swelling, burns, etc.

WARNING

1. People with physical disabilities, sensory or nerve defects, mentally impaired, or young children who do not have a guardian to supervise or guide the use of this massager should not use it.
2. Supervise minors/children when using at all times.
3. Do not allow non-professional repair technicians to disassemble and attempt to repair this massager.
4. Do not use in a humid environment.
5. Do not expose to water or any liquids.
6. Do not place or use the massager near any flammable items.
7. Avoid damaging the massager with sharp objects.
8. Do not expose the massager to extremely high temperatures.
9. Do not stand or walk while the massager is in use.
10. Do not use while sleeping.
11. Do not aggressively pull down on the power cord when unplugging the adapter.
12. Do not plug or unplug the adapter or air pipe connector with wet hands.
13. The massager has a heating surface, make sure to use only as tolerable.
14. If you feel unwell during use or if massager causes discomfort, stop using immediately.
15. It is suggested that each use time should not exceed 15 minutes. Allow the massager to rest 30 minutes before reusing after each session.

IMPORTANT: THE MASSAGER SHOULD ONLY BE USED WHILE SEATED.


PRODUCT SPECIFICATIONS

- Size of product itself: 71*30*28 CM
- Net weight of product itself: 1.87KGS
- Single package size: 34*17.5*25 CM
- Net weight of product including package: 2.1
- Materials: Mesh Fabric + ABS
- Make these points below each bullet points and include in same list as above:

Rated voltage 12V
Rated power 24W

PRODUCT OVERVIEW

1. Foot Region
2. Calf Region
3. Knee Region
4. Thigh Region
5. Adapter
6. Cable Connector Tube
7. Groove
8. Power On/Off
9. Heat
10. Massage Intensity
11. Massage Mode
12. Operation Indication Panel

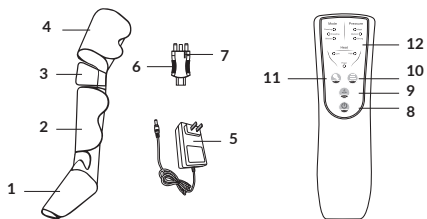

HOW TO USE




1. Put the leg into the correct foot cover. The R on the leg massager represents the right leg and the L on the leg massager represents the left leg.




2. Fit the foot part of the leg sleeve to your foot and strap into place.




3. Fit the calf part of the leg sleeve to your calf and strap into place.

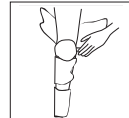


4. Fit the knee part of the leg sleeve to your knee and strap into place.

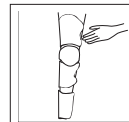


5. Fit the thigh part of the leg sleeve to your thigh and strap into place.

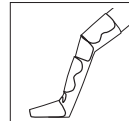


6. Move your legs to make sure both leg sleeves are strapped securely.

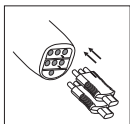


7. Insert the cable connector tube into the corresponding socket of the leg massager.

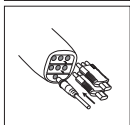

8. Insert the power plug into the corresponding socket of the leg massager.

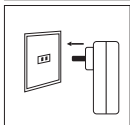

9. Plug the power adapter into a power source outlet.

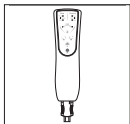

10. Press the power button on the remote control to turn on the leg massager and select your desired massage mode, intensity and heat.

REMOTE CONTROL OPERATION



Pattern	Button	Feature & Operation
	Power	Turn on / off power
	Heat	Heating function: low or high
	Mode	Pressing, kneading or shiatsu
	Intensity	High, medium or low

- Short press the power button to turn on the massager. The default massage setting is pressing mode with low intensity.
- Short press the heat button to turn on the low temperature heating, then press the heat button a second time to turn on high temperature heating. To turn off the heating function, press the heat button a third time.
NOTE: The knee region contains the heating surface.
- Short press the massage mode button. The mode will switch between pressing, kneading, and shiatsu modes of massage.
- Short press the intensity button. The massage intensity pressure will switch from low, to medium and then high.
- Press the power button to stop the massage after use.
- The massager has a set default timer of 15 minutes. It will automatically stop operating after 15 minutes and enter a state.

CLEANING & MAINTENANCE

Please remember to turn off the power before cleaning.

- Clean the massager with a slightly damp towel and then dry it with a soft cloth.
- Do not use detergent, a rough cleaning brush, gasoline, kerosene, glass cleaner or thinner to clean the massager.
- After use and cleaning, store the massager in its original box or in a safe, dry and ventilated place to avoid exposing the massager to moisture or direct sunlight.

TROUBLESHOOTING

If you encounter problems during use, please refer to the following tips and possible solutions.

Problem	Cause	Possible Solution
Does not Operate	The power adapter is not connected to the controller.	Turn on / off power
	The power adapter is not plugged into the socket tightly.	Plug the AC plug of the power adapter into a household outlet tightly.
Automatically turns off	The 15-minute time set by the automatic timer has expired.	Turn on the massager again.
Uneven pressure on the left and right legs	bent tube	Prevent the tube from bending.
Does not heat up	Intensity	Firmly insert the tube connector into the controller.

LIMITED LIFETIME WARRANTY

Your Medic Therapeutics product is backed by a limited lifetime manufacturer's warranty. Medic Therapeutics will repair or replace your device at any time should it fail due to a defect in material or workmanship, subject to the certain limitations.

This limited warranty does not cover any damage that results from unauthorized or improper use, service, or repair. Further, it does not cover damage caused by accident, impact, negligence, or normal wear and tear. Should you discover your Medic Therapeutics product is not functioning properly, please send your device to our repair center for evaluation. If your product cannot be repaired or serviced, we will reserve the right to change it for a similar or newer model.

Please note that a flat rate of \$40.00 will be charged to cover a service evaluation fee and return shipping of your device. All warranty claims must be accompanied by a copy of your proof of purchase from an authorized retailer. Please send your device, proof of purchase, and a check or money order in the amount of \$40.00 made out to Medic Therapeutics to:

Address:
Medic Therapeutics Service Center
3069 Taft Street
Hollywood, FL 33021

Contact:
warranty@medictherapeutics.com