



LED FACIAL MASK

Clinically Tested Light Therapy

USER GUIDE



BE CONFIDENT IN YOUR COSMETICS

L(A)B is the creation of Oxford University Scientists and Sisters Lauren, Charlotte and Sophie with the aim to improve the effectiveness of your everyday skincare and healthcare regimes through cutting edge scientific devices that fit easily into your life.

L(A)B is committed to delivering products of the highest innovation and quality from concept to packaging. We are a unique brand, the first British, Family run Electronic Skincare + Healthcare company. L(A)B are the Beauty Tool Experts.

L(A)B LED Facial Mask

LED Light Therapy Facials tailor made to your exact skincare requirements.

Completely wireless and remote controlled light therapy at home. With a menu of 5 Pre-set facials that are self-timing for 10 or 15 minutes to choose from.

Combining the benefits of 5 different light wavelengths: Red, Blue, Green, Purple and Yellow to tackle all your aging and skincare concerns.

Visibly reduce the appearance of wrinkles, improve skin tone, elasticity and luminosity.


175 specially designed and positioned LED lights cover the entire face and focus on key areas that show signs of aging.

Light weight, and contoured design with adjustable headband is comfortable to wear even when sitting or standing and does not feel heavy on your face.


Soft touch silicon protective eye covers for a comfortable fit and protects eyes so you can even watch TV while wearing your LED Facial Mask.

Product Features

- 1 On/off power button
- 2 Light weight, wireless, contoured LED facial mask
- 3 LED lights
- 4 USB charging port
- 5 Soft touch silicon eye protection
- 6 Adjustable headband
- 7 Chin rest

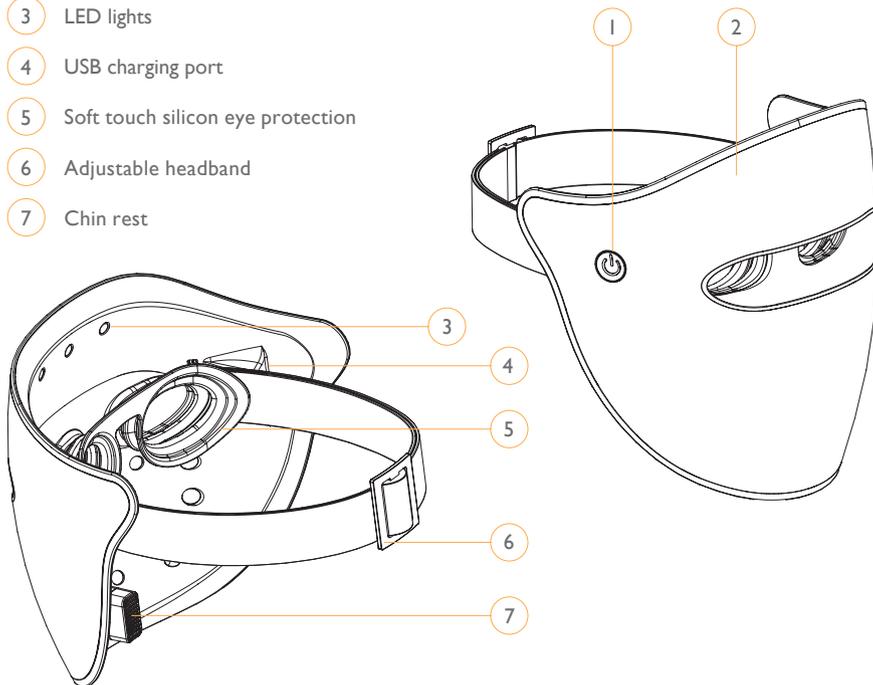


- 8 Digital Screen
- 9 Control Button
- 10 Battery Cover

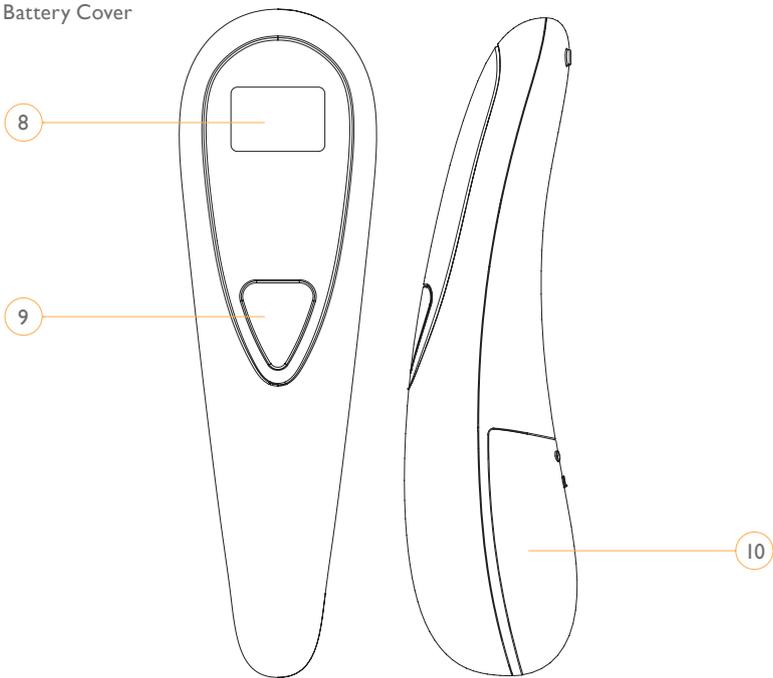

I. Anti-age Facial

Red LED Light + **Skincare**: Visibly reduce the appearance of lines and wrinkles, tighten pores and increase skin firmness, elasticity and bounce back.

1. Start with clean skin.
2. Apply a thin layer of Serum or cream.
3. Put the LED Facial Mask on and adjust the head band so the mask fits comfortably.
4. Press the on/off button on the side of the mask, you will hear a beep and the on/off button will flash green.
5. Press remote control button once to start the Anti-age facial; you will hear a beep and the LED lights will glow red.
6. Relax for 10 minutes. The number of minutes shown on the digital screen of the remote control will count down.
7. When the facial is complete you will hear 2 beeps and the LED lights and remote control will turn off.
8. Remove mask.

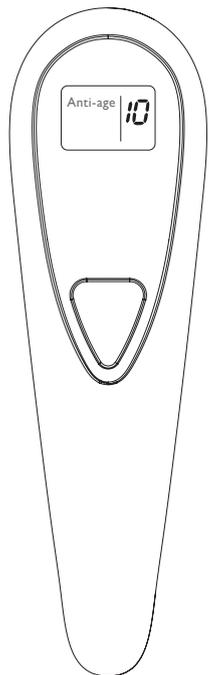


Use 3 times per week.

2. Purify Facial

Blue LED Light: Refines and reduces oiliness and congestion. Helps to re-surface and hydrate and improve skin clarity for a fresh-faced complexion.

1. Start with clean skin.
2. Put the LED Facial Mask on and adjust the head band so the mask fits comfortably.
3. Press the on/off button on the side of the mask, you will hear a beep and the on/off button will flash green.
4. Press remote control button twice to start the Purify facial; you will hear a beep and the LED lights will glow blue.
5. Relax for 10 minutes. The number of minutes shown on the digital screen of the remote control will count down.
6. When the facial is complete you will hear 2 beeps and the LED lights and remote control will turn off.

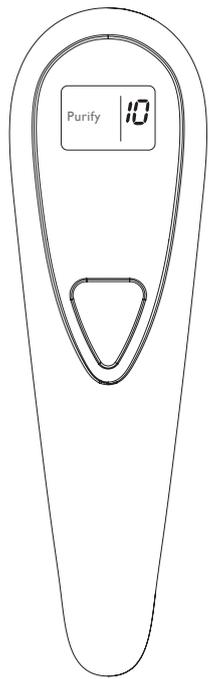


Use 3 times per week.

3. Radiance Facial

Green LED light + **Skincare**: Improve skins' luminosity and radiance. Help to even skin tone and reduce imperfections for a youthful glow.

1. Start with clean skin.
2. Apply a thin layer of **serum or cream**.
3. Put the LED Facial Mask on and adjust the head band so the mask fits comfortably.
4. Press the on/off button on the side of the mask, you will hear a beep and the on/off button will flash green.
5. Press remote control button three times to start the Radiance facial; you will hear a beep and the LED lights will glow green.
6. Relax for 10 minutes. The number of minutes shown on the digital screen of the remote control will count down.
7. When the facial is complete you will hear 2 beeps and the LED lights and remote control will turn off.
8. Remove mask.

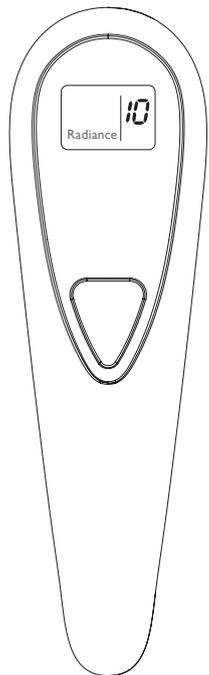


Use 3 times per week.

4. Timeless Facial

Combining the benefits of Red and Green LED lights in a tailor-made pre-set LED facial. To reduce the appearance of wrinkles, hydrate and illuminate the complexion. Help to even skin tone and reduce discolouration and imperfections to give a flawless radiant complexion.

1. Start with clean skin.
2. Apply a thin layer of serum or cream.
3. Put the LED Facial Mask on and adjust the head band so the mask fits comfortably.
4. Press the on/off button on the side of mask, you will hear a beep and the on/off button will flash green.
5. Press the button on the remote control four times until Anti-age and Radiance are visible. The LED lights will glow yellow and the pre-set facial will start.
6. Yellow Light = Anti-age + Radiance for 5 minutes.
Then, Red Light (Anti-age) for 5 minutes.
Then, Green Light (Radiance) for 5 minutes.
7. When the treatment is complete you will hear 2 beeps and the LED lights and remote control will turn off.
8. Remove mask.

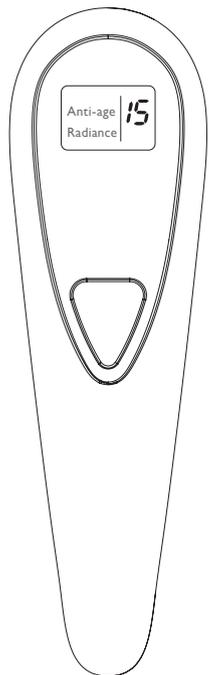


Use 3 times per week.

5. Balance Facial


Combining the benefits of Red and Blue LED lights in a tailor-made pre-set LED facial. To reduce the appearance of wrinkles and the signs of aging, and balance the complexion reducing oiliness and dryness. Refining pores and clarity to give an even, balanced and youthful complexion.

1. Start with clean skin.
2. Apply a thin layer of serum or cream.
3. Put the LED Facial Mask on and adjust the head band so the mask fits comfortably.
4. Press the on/off button on the side of mask, you will hear a beep and the on/off button will flash green.
5. Press remote control button five times until Anti-age and Purify are visible. The LED lights will glow purple and the pre-set facial will start.
6. Purple Light = Anti-age + Purify for 5 minutes.
Then, Red Light (Anti-age) for 5 minutes.
Then, Blue Light (Purify) for 5 minutes.
7. When the facial is complete you will hear 2 beeps and the LED lights and remote control will turn off.
8. Remove mask.

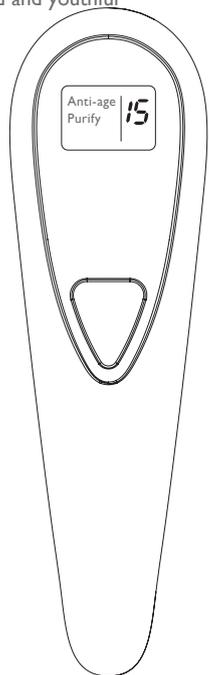

Use 3 times per week.

Charging

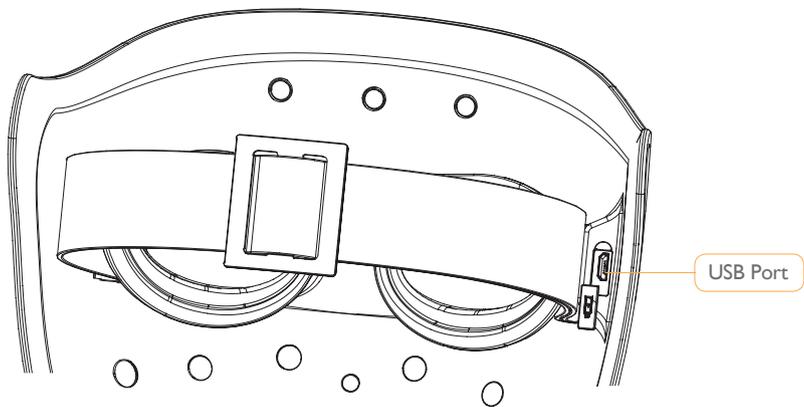

1. Connect the micro-USB in to the charging port on the mask (see above).
2. Connect the USB end of the charging cable in to device/ adaptor plug.
3. On/off button will flash green while charging.
4. Full charge takes approximately 6 hours.
5. When fully charged the on/off button light will glow steady green.
6. Full charge will last approximately 100 minutes or 6-10 facial treatments.
7. When the battery is low and requires charging, the on/ off button will flash red.
8. Remote control battery will last approximately 600 minutes before needing replacement.

Frequently Asked Questions

1. Will the Lights get hot?

No, they are cool LED lights, that do not heat up.

2. Are the lights safe for my eyes?

Yes, they are. They are LED lights, not lasers.

3. Can I do more than one facial after the other?

Yes.

4. Is my mask turned on?

The on/off button will have a green light around it.

5. How often do I need to charge my mask?

Depending on usage the mask will need charging every 6-10 facials. When the mask is low on battery and requires charging the on/off button will flash red.

6. Is the mask heavy?

The mask is designed to sit lightly on the face so you can wear while sitting up and even watching TV.

7. Can I stop a facial early?

Yes, if you want to stop a facial you can remove the mask, press the on/off button on the side and press the remote control button until the screen goes blank.

8. Can I switch between facials?

Yes, when you change between facials the timer will start again on the remote control screen.

You can also carry out different facials throughout the week or at different times of the year to tailor-make them to your requirements.

Quick Start Guide

1. When you receive mask put on charge until the on/off button glows a steady green.
2. Open battery cover on remote control, remove the plastic disc from the battery, and close the battery cover leaving the battery in place.
3. Try mask on and adjust the head band until the mask fits comfortably and securely on your face.
4. When ready to start an LED facial, start with clean skin and choose which facial you would like to use.
5. Prepare skin with relevant skincare product, then put on the LED face mask.
6. Press on/off button on side of mask, you will hear a beep.
7. Press the remote control button until the required facial has been selected, the name will be shown on the screen.
8. Relax (the minutes will count down on the screen of the remote control).
9. When the facial is finished you will hear 2 beeps and the LED mask and remote control will turn off.
10. Remove the LED Facial mask, your treatment is complete

CAUTIONS: Do not use the LED Facial mask if you suffer from light sensitive medical conditions including photo-sensitive epilepsy, if you have sun burnt or inflamed skin, if you are pregnant. Do not stare directly at the LED lights.

To clean wipe the mask with a clean, damp, tissue or cloth and towel dry. Do not submerge in water.



Contact us

E: info@lab-lifeandbeauty.com | T: 01722 710 811 (option 3) UK customer services

www.lab-lifeandbeauty.com

