


NAMES AND FUNCTIONS OF THE CONTROL



POWER BUTTON:

Press to turn the product on or off.

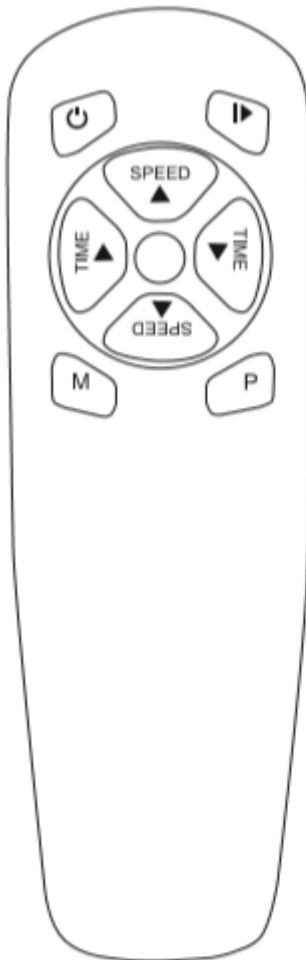

START/STOP:

Plug in the product, use the power button to turn on and off, then use this button to start and stop the product.

TIME:

The default setting is 10 minutes. Press the Time button to increase or decrease the time.

SPEED:

Press this button to increase or decrease the speed level.



M:

This is the Speed Quick Mode button for levels 3, 6 and 9. To access each speed level, press once for level 3, twice repeatedly for level 6, and 3 times repeatedly for level 9.



AUTO MODE:


Press this button to access the automated setting. Note: Time nor Speed can be adjusted under this setting.




Reception range for the controller is 8 feet. Ensure the infrared red emission indicator on the controller is aligned with the receiving window on the product.

COMPONENTS AND USE

1. Plug in the AC adapter.
2. Tightly connect the main body of the platform with the plug on the side of the platform.
3. Place your foot or legs onto the massager to relieve pain, tensions and fatigue. Use the foot massager to your appropriate needs or liking.
4. Unplug when finished using.

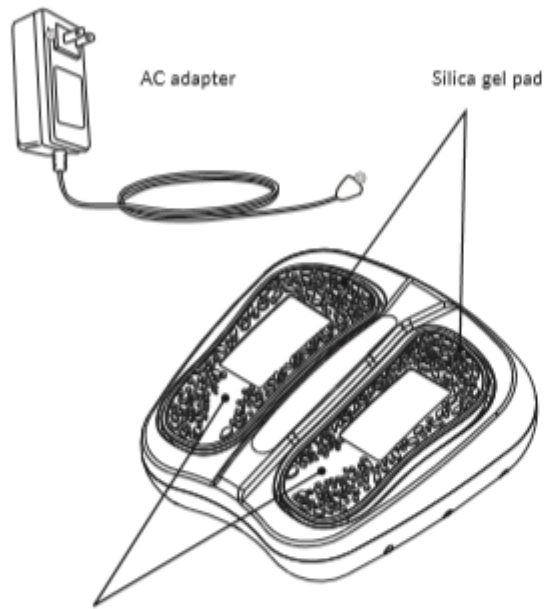

Foot massage structure. Place shank or arch of the foot on the machine (Take off Silican gel pad, can be cleared, thoroughly dry to use after cleaning.)

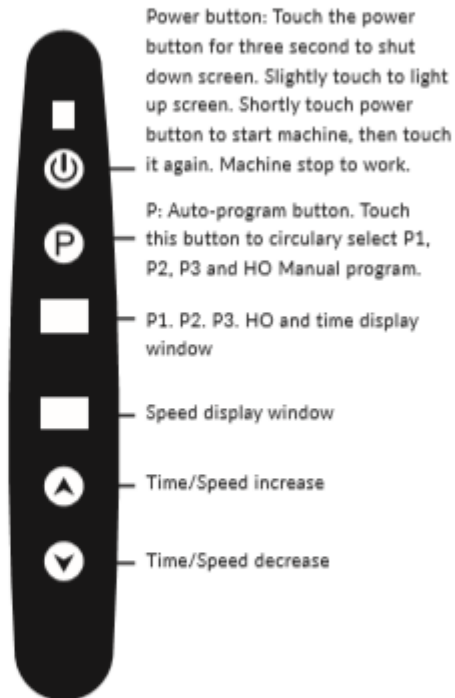

CORRECT USAGE

