




MEDIC THERAPEUTICS MESSAGE GUN

Operating Instructions

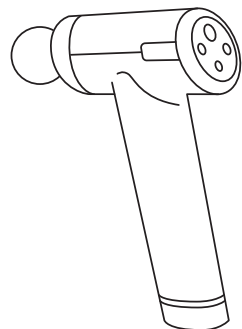

2400r/min

Multi-range Adjustment


Power Display

Deep Amplitude 12MM

Switch Control

Durable Life

Please read carefully before using

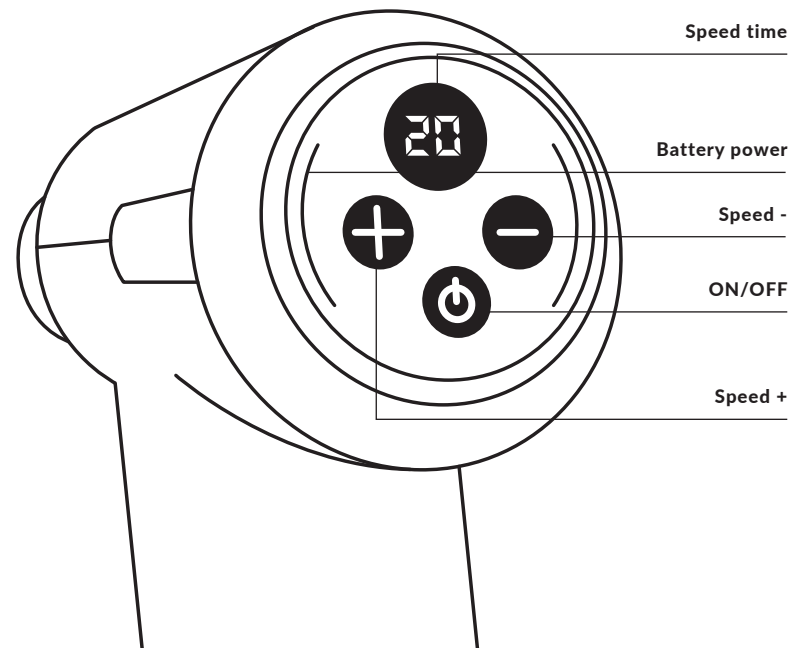

1. SPECIFICATIONS



U-Shaped Head
(For use on either side of the spine)



Shovel Head
(For full body use)



Bullet Head
(For smaller, localized use)




Large Round Head
(For larger muscle areas)



Flat Head
(For use on legs, hips and hands)



Mushroom Head
(For use on tender areas)

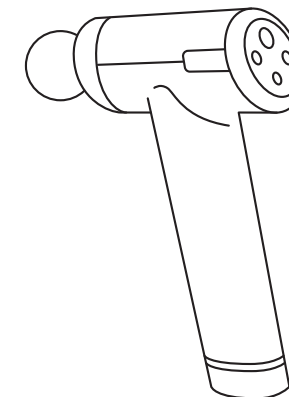

Vibration Frequency:
20Levels LEVEL1-20:30-50HZ
Product Weight: 0.70kg
Product Size: 175*65*238mm
Charge: 110-240v 50/60Hz DC,12.6V=1A
Main Material: Nylon (Body of Massage gun, Massage head)
Circle-shaped Massage head's material: EVA lithium battery
Operational environment: -10~40°C

2. INSTRUCTIONS FOR USE

Charging:

- Prior to use, charge the device for up to 6 hours.
- To charge, connect the DC end of the adapter (included) to the battery's charging port and plug the AC end into the wall outlet.
- The LED battery band light indicates remaining battery life and pulsates when device is actively charging.
- Device is fully charged when the LED light stops pulsing.
- Device may be recharged at any battery level.
- It is not recommended to allow the battery to lose all power before recharging.
- Depending on speed level and pressure, the battery should last 2+ hours.

Operating:

- Ensure the device is locked into place in the provided base.
- To turn on the device, move the switch at the bottom of the battery to the "ON" position. Move switch to the "Off" position to turn the device off.
- To adjust settings, once the device is switched on, press the power button once for Level One twice for Level Two, and again to switch the device off. The power level will be indicated by the blue LED light.
- Massage body, applying pressure when desired.
- To remove the battery, lightly press the battery release button. Then firmly pull downwards on the battery.

Maintenance, Cleaning and Storage:

- Clean with a slightly damp cloth. Dry with a soft cloth.
- When traveling, remove the battery and store in the supplied carton box.

Warnings:

- Keep away from liquids or heat sources and keep ventilation ports free from dust and debris.
- Do not immerse in water.
- Do not remove screws or attempt to disassemble.
- Do not operate for more than one hour at a time. Allow 30 minutes in-between uses.
- Unplug the charger from the battery after charging and prior to use.
- Do not continue use if injured and consult your physician.

Safety Instructions:

For adult use only.

- Use only on clean, dry areas of the body, over clothing. Lightly press, moving across the body for approx. 60 seconds per region.
- Use the Hypervolt function only on soft tissue. Do not use on the head or any hard or bony part of the body.
- Use only the applicators that provide the desired effect.
- Bruising may occur regardless of setting or pressure. Frequently check treated areas and discontinue use at the first sign of pain or discomfort.
- To avoid pinching, keep hair, fingers and any other body part away from the shaft and backside of the device.
- Do not place any objects into the ventilation holes.
- Only recharge with the supplied charger.
- Carefully examine the device and battery before each use.
- Do not tamper with or alter the device in any way.
- Never leave the device unattended while switched on or charging.

3. DISCLAIMER

If any of these medical conditions apply to you, consult with your physician before using the Massage Gun: Pregnancy, Pacemaker, Diabetes with complications, IUD, recent joint replacement, metal pins or plates, herniated disks, or, spondylolysis. Frail individuals and children should be accompanied by an adult with using the device.

Ongoing research for vibrational therapy and exercise related to specific medical conditions indicates that the above list will soon be shortened. Practical experience has shown that vibrational therapies and exercise should be incorporated in many treatment plans. This must be done on the advice of a physician, specialist, or physiotherapist.

4. NON MEDICAL USE

Medic Therapeutics is a trademark name for a family of products that includes a variety of emergency and wellness supplies such as massage guns, thermometers, face coverings and masks (such as coverings of the nose and mouth that may include masks as well as bandanas, scarves, and other cloth coverings). None of Medic Therapeutics products are designed or intended to prevent, diagnose or treat any illness, condition or disease, even when used as directed. This product is not a medical device. This product is not intended for use in the diagnosis of disease or other conditions or in the cure, mitigation, treatment, or prevention of disease. Always consult your doctor or physician for any medical questions you may have.

This massage gun is designed to relieve body stress and pain by manipulating soft tissues in the body. Always consult your doctor or physician immediately for any medical questions you may have. This apparatus is NOT guaranteed as a prevention or treatment for any disease or medical condition. For proper use only. Operate only as directed. Misuse may result in sickness or death. Always carefully review instructions for use in the operating instructions manual prior to use.

Keep out of reach of children and supervise children under the age of eighteen (18) at all times if the minor is utilizing the massage gun. Do not allow the user to place this item in his/her mouth at any time. No child or minor should ever be left unattended with any of these items, as small parts may become loose & pose a choking hazard.

5. LIMITED LIFETIME WARRANTY

Your Medic Therapeutics Massage Gun is backed by a limited lifetime manufacturer's warranty. Medic Therapeutics will repair or replace your device at any time should it fail due to a defect in material or workmanship, subject to the certain limitations.

This limited warranty does not cover any damage that results from unauthorized or improper use, service, or repair. Further, it does not cover damage caused by accident, impact, negligence, or normal wear and tear. Should you discover your Medic Therapeutics Massage Gun is not functioning properly, please send your device to our repair center for evaluation. If your product cannot be repaired or serviced we will reserve the right to change it for a similar or newer model.

Please note that a flat rate of \$35.00 will be charged to cover a service evaluation fee and return shipping of your device. All warranty claims must be accompanied by a copy of your proof of purchase from an authorized retailer. Please send your device, proof of purchase, and a check or money order in the amount of \$35.00 made out to Medic Therapeutics to:

Medic Therapeutics Service Center
3069 Taft Street
Hollywood, FL 33009



*This Prop 65 warning label is affixed to this product as mandated by the State of California relating to all consumer products. **WARNING:** This product may contain chemical(s) determined by the State of California to cause of cancer, birth defects and/or reproductive harm. For more information go to <https://www.p65warnings.ca.gov/>*