



HEATED VEST

USER MANUAL

SAFETY PRECAUTIONS:

VEST:

1. Read user manual prior to first use.
2. It is not recommended to allow the heating elements of the vest to come into contact with bare skin.
3. Never operate the vest with heating settings if the inner liner is wet.
4. Shut off the power immediately if any discomfort occurs during use.

POWER BANK:

1. Close supervision is necessary when the power bank is used near children to reduce the risk of injury.
2. Do not place your fingers or hand inside the product.
3. Do not expose the power bank to rain or snow.
4. Use of a charger or power supply charger not sold or recommended by the power bank manufacturer may result in the risk of fire or injury.
5. Do not use the power bank in excess of its output rate. Overloading outputs above the total rate may result in the risk of fire or injury.
6. Do not disassemble the power bank. When service or repair is required, take the battery bank to a qualified service representative. Incorrect reassembly may result in the risk of a fire or injury.
7. Have servicing performed by a qualified representative using only identical replacement parts. This will ensure that the safety of the product is maintained.

PRODUCT SPECIFICATIONS:

VEST:

- Vest Material: Shell 100% polyester, Lining 100% polyester
- Battery: Power Bank operated (Output: 5V 2.1A)
- Heating: Carbon fiber heating elements
- Temperature Settings: High, Medium & Low

POWER BANK:

- Material: ABS
- Unit Size: 146*71*16.1mm
- Battery Capacity: 10000mah Li-polymer battery
- Battery Life: Up to 500 cycle times or more
- DC Input: 5V 2A(max)
- DC Output: 5V 2A(max)
- Accessory: Micro USB Charging Cable, Manual



OPERATING INSTRUCTIONS:

1. Put on the vest and grab the USB charging lead found in the left inner pocket. Plug the USB lead into the provided power bank. Turn on the power bank and then insert the bank into the pocket. (Power Bank- Output: USB 5V 2A, Input: Micro 5V 2A)
- Charging the Power Bank: Connect the USB port of the supplied micro charging cable to an AC powered USB adaptor, computer or laptop. Then, connect the micro USB port of the same cable to this power bank. Once connected, the 4 white lights will flash, indicating that charging has begun. Once fully charged, all four lights will stay illuminated consistently.
- Use Power Bank to Charge Your Devices: Connect the USB Port of your phone charging cable to the USB output port of power bank. Connect the other end of your charging cable to your phone. Once connected to your phone, charging commences immediately. The 4 white power indicator lights will blink off accordingly as the battery power is consumed.

- Press and hold the power button of the vest for 3-5 seconds to turn the power on/off. Heating will commence.
- Press and release the button to alternate between the 3 heating settings. The light of the power button will illuminate a different color, indicating whether the vest is set to High, Medium or Low Heat. Allow 25-30 minutes for the vest to warm up to the desired heating setting.

a. Red: (High)



b. White: (Medium)



c. Blue: (Low)




Note: When first turning on the power button, the power button will illuminate red; it will take 5 minutes for the vest to warm up to the designated high heat temperature. After the first 5 minutes have passed- if you have not switched between heating settings during this time -the vest will automatically initiate Overheating Protection and switch to the medium heat (the power button will illuminate white).


- To stop the heating, press and hold the power button until the vest shuts off, or unplug the USB charging lead from the power bank.

HEATING AREAS:

Note: When the power bank is fully charged, vest heating lasts for approximately 6 hours (when set to LOW heat).

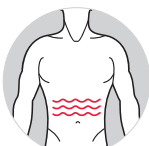

Abdomen

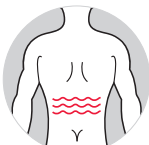

Back



Waist

WASH & CARE INSTRUCTIONS:

- Hand wash only
- Water temperatures should not exceed 86°F
- Do not bleach. Do not tumble dry
- Do not dry clean
- Do not iron
- Remove power bank prior to washing

LIMITED LIFETIME WARRANTY:

Your Medic Therapeutics Heated Vest is backed by a limited lifetime manufacturer's warranty. Medic Therapeutics will repair or replace your device at any time should it fail due to a defect in material or workmanship, subject to the certain limitations.

This limited warranty does not cover any damage that results from unauthorized or improper use, service, or repair. Further, it does not cover damage caused by accident, impact, negligence, or normal wear and tear. Should you discover your Medic Therapeutics Heated Vest is not functioning properly, please send your device to our repair center for evaluation. If your product cannot be repaired or serviced, we will reserve the right to change it for a similar or newer model.

Please note that a flat rate of \$35.00 will be charged to cover a service evaluation fee and return shipping of your device.

All warranty claims must be accompanied by a copy of your proof of purchase from an authorized retailer. Please send your device, proof of purchase, and a check or money order in the amount of \$35.00 made out to Medic Therapeutics to:

Address:
Medic Therapeutics Service Center
3069 Taft Street
Hollywood, FL 33021

Contact:
warranty@medictherapeutics.com