


SUGGESTED TREATMENT AREAS



Wrist



Knee

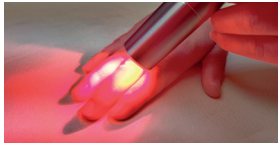

Fingers

LIMITED LIFETIME WARRANTY

Your Medic Therapeutics product is backed by a limited lifetime manufacturer's warranty. Medic Therapeutics will repair or replace your device at any time should it fail due to a defect in material or workmanship, subject to the certain limitations.

This limited warranty does not cover any damage that results from unauthorized or improper use, service, or repair. Further, it does not cover damage caused by accident, impact, negligence, or normal wear and tear. Should you discover your Medic Therapeutics product is not functioning properly, please send your device to our repair center for evaluation. If your product cannot be repaired or serviced, we will reserve the right to change it for a similar or newer model.

Please note that a flat rate of \$38.00 will be charged to cover a service evaluation fee and return shipping of your device. All warranty claims must be accompanied by a copy of your proof of purchase from an authorized retailer. Please send your device, proof of purchase, and a check or money order in the amount of \$38.00 made out to Medic Therapeutics to:

Address:
Medic Therapeutics Service Center
3069 Taft Street
Hollywood, FL 33021

Contact:
warranty@medictherapeutics.com



RED LIGHT THERAPY



User Manual

Please read this manual prior to use in order to guarantee the safe and proper use of this device.

ABOUT RED LIGHT THERAPY

Red light therapy is intended to assist in pain relief and alleviation of inflammation. Our red light therapy device uses premium tri-spectrum technology to help these symptoms recover through active cell regeneration; our device uses a combo of red 630nm, 660nm and 850nm infrared LEDs which are believed to have a deep penetration effect.

Red light therapy is believed to have many benefits such as:

- Anti-aging & anti-inflammatory effects on the skin
- Stimulating fibroblasts to produce collagen and elastin
- Providing essential energy for cellular renewal
- Improving skin strength, elasticity and plumpness and reducing arance of fine lines and wrinkles
- Increasing circulation
- Assisting with tightening and firming of the skin
- Assisting with lymphatic drainage to help eliminate toxins
- Stimulating production of cellular ATP
- Temporarily relieve arthritic pain, muscle and joint stiffness, inflammation and swelling

WARNING

1. Do not use this Device on infants and young children.
2. Do not shine Device lights directly onto eyes.
3. Discontinue use if Device overheats or becomes uncomfortably hot to your skin.
4. Do not use Device if your skin is light sensitive or you are currently using medication that may cause sensitivity of skin to light.
5. Ask your Doctor prior to using if you:
 - Are pregnant or suspect you may be pregnant
 - Are taking drugs known to cause light sensitivity
 - Have any suspicious or cancerous lesions
 - Hemorrhagic diseases
6. Contact your Doctor if you have any medical questions regarding use of Device. Contact your Doctor if you suspect that you have had a negative reaction from using this Device.

PRODUCT SPECIFICATIONS

- Size of product itself: D26*120.5mm
- Net weight of product itself: 0.21KG
- Single package size: 163(L)*53(W)*226(H)mm
- Net weight of product including package: 0.5g
- Material: Stainless-Steel
- Battery (QTY & Type): 2pcs & 18650


INCLUDED IN PACKAGING

- 1x Red Light Therapy Device
- 1x Wrist band
- 2x 18650 rechargeable batteries
- 1x USB Cable
- 1x Protective Eyewear
- 1x Carrying bag

HOW TO USE

Battery Installations:

- Insert the battery provided into the battery compartment: positive pole against the TOP and negative pole against the BOTTOM.
- Only use the included 18650 battery.

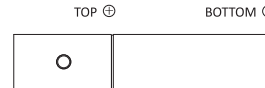

HOW TO OPERATE

- Press the button located on the top of the device to turn it on. The red light will illuminate, and the therapy effects will be activated. The default operating time is 5 minutes (after 5 minutes, the device will automatically shut off).
- Hold the unit over the painful area at a minimum distance of 1 inch away from the skin. Keep the unit hovering over the treatment area. The red light therapy device penetrates up to 1 inch into your tissues/joints.
- Hold the unit over the treating area for 5 minutes, guides as below:

Pain Level	Light Pain	Average Pain	Deep Pain
Treatment	Treat Spot Once a Day	Treat Spot Twice or 3 Times a Day	Treatment Spot 3 Times a Day

- Do not keep the device on hovering over an area for more than 5 minutes at a time. If you want to treat several parts of your body consecutively, TURN OFF the device for at least 1 minute after every 5 minutes of operation. This will keep the unit from overheating.

CHARGING THE DEVICE

Please note, while charging, the charging indicator light of the device will be red. Once fully charged, this indicator light will switch to green.