



**MEDIC
THERAPEUTICS**

CORDLESS HAND MASSAGER WITH HEAT



User Manual

SAFETY PRECAUTIONS



Please read the following instructions carefully before putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these operating instructions for further use!

- The massage device conforms to recognized technical principles and the latest safety regulations.
- Do not wet, do not use pins, never remove cover.
- This item is NOT A TOY. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
- This appliance should never be left unattended when plugged in.
- Any possible repairs may only be carried out by authorized expert staff. Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- Do not use any damaged cables, plugs or loose sockets.
- Never operate under blankets where the air opening may be blocked.
- If plugs or cords are damaged, they must be replaced by the manufacturer, a service representative or by qualified personnel.
- In case of malfunction, disconnect immediately from the mains.
- Do not use if you have skin disorders, open wounds, or swollen or inflamed areas.
- Misuse or incorrect use excludes any liability for damage.
- Do not use while you sleep.
- Do not use this product for medical treatment.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time should not be exceeded 15 minutes at a time.
- The packaging material supplied should not be used as a toy.
- Consult a doctor prior to use if you are one of:
 - Those who are using electronic medical instrument imbedded in the body, such as cardiac, pacemaker etc.
 - Those who are being treated by doctors, especially those who feel discomfort.
 - Those who are malignant tumor patients, heart disease patients, or acute patients.
 - Women who are pregnant or menstruating.
 - Those with osteoporosis or a fractured spine.
 - Those with skin disease or those whose skin is impaired.
- Failure to observe the above instructions may constitute misuse of the product and can cause serious injury or burns.

PRODUCT FEATURES

- Compression & Vibration Massage Functions
- 3 Different Intensity Levels
- Rechargeable battery included
- 15 Minute automatic timer
- Lightweight & easily portable

PRODUCT OVERVIEW

- Hand Massager with Heat
- 1x USB Cable



OPERATION INSTRUCTIONS & CONTROL PANEL OVERVIEW

ON/OFF Button: Press and hold for 2 seconds to turn on the device. The ON/OFF button light will turn green and activate the device to work on default Mode 1 (Green Indicator Light) at Medium intensity level (Yellow Indicator Light). To turn off, press and hold the power button for 2 seconds.

Pressure Mode button:

Press to set the Pressure Mode at Mode 1, Mode 2 or Mode 3

Indicator:


Mode 1- Green Light
Mode 2- Yellow light
Mode 3 - Red light

Vibration Mode Button:

Press to activate vibration function. You may cycle through Mode 1, Mode 2 or Mode 3

Indicator:

Mode 1- Green Light
Mode 2- Yellow light
Mode 3 - Red light

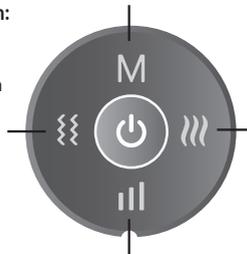

Heat Button:

Press to activate heat function. You may cycle through Low, Medium, High or Off settings.

Indicator:

Low - Green light
Medium - Yellow light
High - Red light

Intensity Button:

Press to set the massage intensity at Low, Medium, or High Level.

Indicator:

Low - Green Light
Medium - Yellow light
High - Red light

The Indicator light flashes when the battery is low.

CHARGING INSTRUCTIONS

1. Connect the USB cable to the corresponding port on the device.
2. Normally it takes 2.5-3 hours to fully charge the battery.
3. While charging, the indicator light turns red, when fully charged, the indicator light will turn green.
4. If the battery is full, the device can work continuously for around 1.5 hours.
5. You are able to use the massager while charging, but it is recommended to charge the battery first then use.



- This device will switch off automatically after 15 minutes.
- Do not setup or use the device in a bathroom or similarly wet / damp area.
- Do not use this device while driving.
- You should allow your body some rest periods. To avoid straining your muscles too much, we recommend to not exceed a continuous massage of more than 15 minutes.
- Ensure the device is not exposed to direct sunlight or placed near sources of heat.

CLEANING & MAINTENANCE

1. Dampen a cloth in water or a 3%-5% mild detergent solution.
2. Wipe the dirty areas with the wet cloth.
3. Wait for the unit to become fully dry before use.
4. Always unplug the unit before cleaning it.
5. Before cleaning, allow unit to cool.
6. Wipe the unit with a soft, dry cloth. Do not use cloths containing any kind of chemical or alcohol and solvent liquids.
7. Never submerge any part of the unit in liquid.
8. Place massager in a safe, dry, and cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the surface.
9. Keep away from children at all times.

TROUBLESHOOTING

| Trouble | Reason/Solution |
|--------------------------------|---|
| Cannot start the massager | Battery Low. Please charge it before use |
| The massager stoped abnormally | Battery runs out of power. Charge the battery the Self-timer 15 mins is out,restart the product |
| Massager makes slight noise | It is normal working noise of the internal mechanisms |

PRODUCT SPECIFICATIONS

- Size of product itself: 11 X 7 X 5.5 IN
- Net weight of product itself: 0.75 KGS
- Single package size: 30 X 20 X 14.3 CM
- Net weight of product including package: 1.17 KGS
- Materials: ABS
- Battery Voltage: 3.7VDC 2200mAh
- Fully Charged Battery Time: ≤ 3 hours
- Battery Max. Runtime: ≥ 1.5hours
- Automatic Timer: 15 Minutes

LIMITED LIFETIME WARRANTY

Your Medic Therapeutics Product is backed by a limited lifetime manufacturer's warranty. Medic Therapeutics will repair or replace your device at any time should it fail due to a defect in material or workmanship, subject to the certain limitations.

This limited warranty does not cover any damage that results from unauthorized or improper use, service, or repair. Further, it does not cover damage caused by accident, impact, negligence, or normal wear and tear. Should you discover your Medic Therapeutics Product is not functioning properly, please send your device to our repair center for evaluation. If your product cannot be repaired or serviced, we will reserve the right to change it for a similar or newer model.

Please note that a flat rate of \$30.00 will be charged to cover a service evaluation fee and return shipping of your device. All warranty claims must be accompanied by a copy of your proof of purchase from an authorized retailer. Please send your device, proof of purchase, and a check or money order in the amount of \$30.00 made out to Medic Therapeutics to:

Address:

Medic Therapeutics Service Center
3069 Taft Street
Hollywood, FL 33021

Contact:

warranty@medictherapeutics.com





**MEDIC
THERAPEUTICS**