

Thermapulse™
Relief Wrap

Extra-Long Massaging Heat Wrap

**IMPORTANT SAFETY INSTRUCTIONS
READ BEFORE FIRST USE**



WARNING: THIS PRODUCT CONTAINS MAGNETS. CAN BE HARMFUL TO THOSE WHO HAVE A PACEMAKER, DEFIBRILLATOR OR OTHER MEDICAL IMPLANTS. FAILURE TO FOLLOW THIS WARNING CAN RESULT IN SERIOUS INJURY.

WARNING: This appliance has a polarized plug (one blade is wider than the other). To reduce the risk of electric shock, this plug is intended to fit into a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to modify the plug in any way.

CAUTION: HANDLE WITH CARE. ALWAYS LAY FLAT OR FOLD. CRUSHING OR CRIMPING OF THE HEATING WIRES CAN LEAD TO PREMATURE FAILURE.

IMPORTANT SAFETY INSTRUCTIONS

DO NOT DESTROY


DANGER: TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK AND FIRE, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS.

- 1. READ ALL INSTRUCTIONS.**
- 2. DO NOT USE WHILE SLEEPING.**
- 3. DO NOT USE ON AN INFANT.**
- 4. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, SLEEPING OR UNCONSCIOUS PERSON, OR A PERSON WITH POOR BLOOD CIRCULATION UNLESS CAREFULLY ATTENDED.**
- 5. DO NOT USE ON AREAS OF INSENSITIVE SKIN.**
- 6. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING, CHECK SKIN UNDER PAD FREQUENTLY.**
- 7. NEVER USE PAD WITHOUT COVER IN PLACE.**
- 8. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.**
- 9. DO NOT SIT ON OR CRUSH PAD - AVOID SHARP FOLDS.**
- 10. NEVER PULL THIS PAD BY THE SUPPLY CORD.**
- 11. DO NOT USE THE CORD AS A HANDLE.**
- 12. CAREFULLY EXAMINE INNER COVER BEFORE EACH USE. DISCARD THE PAD IF INNER COVERING SHOWS ANY SIGN OF DETERIORATION.**
- 13. LOOP CORD LOOSELY WHEN STORING, TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS.**
- 14. SAVE THESE INSTRUCTIONS.**

This heating pad has a polarized plug (one blade is wider than the other) as a safety feature. This plug will fit into a polarized outlet only one way. If the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.


USAGE INSTRUCTIONS

STEP 1:

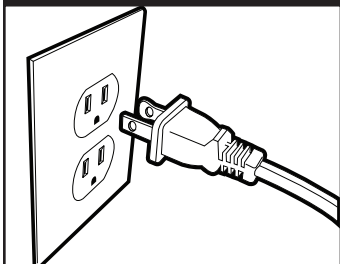


Plug the Relief Wrap™ into a 120V AC electrical outlet.

STEP 2:

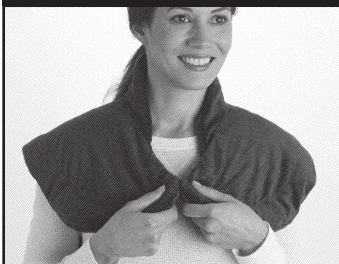

Place the Relief Wrap™ on body as shown so the magnets will connect together.

STEP 3:

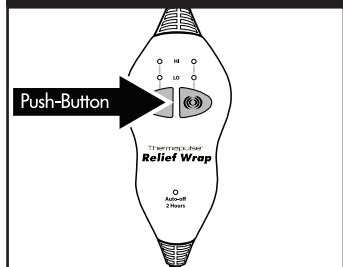

To turn on the Relief Wrap™, select the desired heat or massage function.


IMPORTANT: If the Heat Indicators DO NOT light up or stay on, unplug from electrical outlet.

PRE-HEAT:

Allow several minutes for the Relief Wrap™ to preheat.

IMPORTANT: If there is a red blinking "HIGH" light, this means an error has occurred. Unplug Wrap from outlet, and inspect the condition of the cord and the function of the AC outlet. If blinking continues, stop usage of Wrap and contact Customer Service.


HEAT:

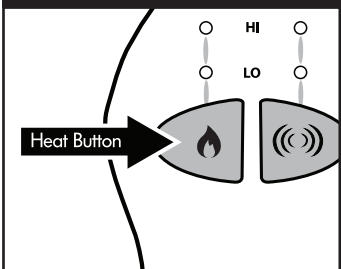

To apply Heat, press the Heat Activation Button and select your desired temperature: Hi or Lo.

Press **1X** for Lo Heat
Press **2X** for Hi Heat
Press again for Off

MASSAGE:

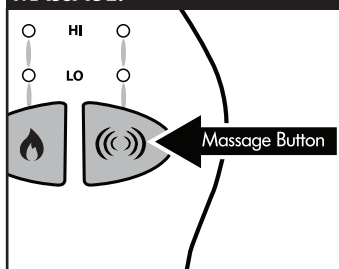


To activate the Massage feature, press the Massage Activation Button and select your desired speed:

Press **1X** for Lo Steady
Press **2X** for Lo Pulse
Press **3X** for Hi Steady
Press **4X** for Hi Pulse

Note: Your Relief Wrap™ is designed for a strong massage. If a lighter massage is desired, simply flip the product and wear it with the controller on the LEFT side. The massagers will feel lighter when worn this way.

AUTO-OFF FEATURE

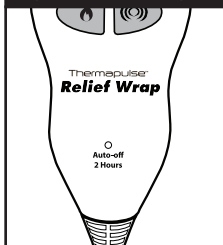

The Auto-Off Feature is included in the Control on the Relief Wrap™. The Auto-Off allows for the Relief Wrap™ to conserve energy by turning off after 2 continuous hours of use. The Massage feature will automatically shut off after 15 minutes of use.

IMPORTANT: The Auto-off timer will reset after a power outage.

WASHING INSTRUCTIONS:

**NOTE: ONLY COVER IS WASHABLE! DO NOT WASH HEATING PAD!
DO NOT DRY CLEAN!**

DO NOT USE BLEACH. DO NOT WRING.

1. Unzip the cover on the side and unfasten the inner straps which connect the heating pad to the cover. (See Fig 7) Then carefully remove cover from Relief Wrap™.
2. Machine wash cover in warm water on gentle cycle, tumble dry on medium heat, remove quickly from dryer.
3. Relief Wrap™ can be hand washed in warm, soapy water. Hang to dry.
4. When dry, place the heating pad back inside the cover and re-fasten the straps zip closed when done.



CAUTION: HANDLE WITH CARE. ALWAYS LAY FLAT OR FOLD. CRUSHING OR CRIMPING OF THE HEATING WIRES CAN LEAD TO PREMATURE FAILURE.

IMPORTANT SAFETY INSTRUCTIONS READ BEFORE FIRST USE

100% Polyester • 26" x 33" • 70 WATTS • 120 VAC ONLY • 60 HZ

FOR HOUSEHOLD USE ONLY



Model #17001

Distributed By



©2017 Ontel Products Corp. • Fairfield, NJ 07004
All Rights Reserved • Made In China
RWIM30517-NKS