



Give your family a warm and healthy life.


Product Features



Products for legs and knees

Product Introduction

These products are suitable for legs, knees and feet.

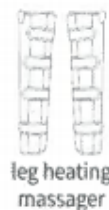


Give your family a warm and healthy life.

Usage Methods

1. Whilst sitting down, legs at 90° angle, place the wraps over each leg and adjust the velcro straps to the comfort of your liking
2. **Double Switches Hand Controller:** Push upwards the switches to turn on heating and vibration. From top to bottom the 4 gears are: High(mode 1), Medium(mode 2), Low(mode 3) and Stop. Please note: the device will auto shutoff after 30 minutes of use.
3. Medium heating (55°C) is recommended, 30-60 minutes for each use, 1-2 times per day. Consuming a glass of water is suggested after each use.

4. If you choose to use the product continually for more than 2 hours, low heat setting is recommended. The High heat setting (65° C) is not recommended for extended use (i.e. for more than 2 hours at a time.) This avoids frequently overheating and extends the product lifetime.

CAUTION: During use, consumer must be aware of heat setting to avoid possible over heating.

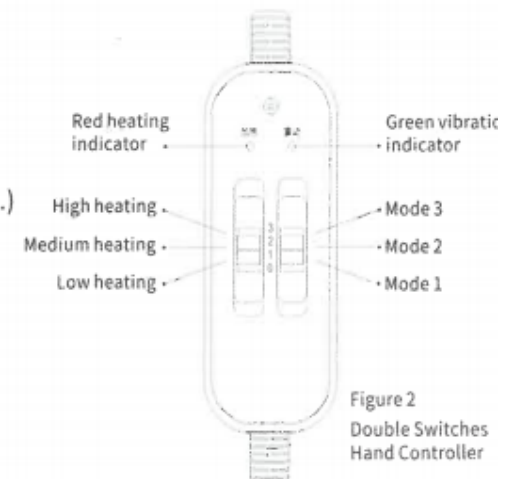



Matter of Attention

1. The product is not waterproof, and for indoor use only. Do not use when wet.
2. The product should not be under heavy pressure and can not be folded for use.
3. Discontinue use if damage to heater or power cord is found.
4. Oil or other flammable substances sticking to heater should be cleaned before use.
5. Do not cut, damage, refit the heater and do not insert sharp objects.
6. Please turn off the heater and disconnect power after use.
7. Do not use during sleep.
8. Disconnect power immediately and stop using it once you feel uncomfortable or unexpected situation occurs during use.
9. Gravidas, diabetics, and malignant tumor patients should consult doctors use of device

Contraindications

1. Patients with severe cardiovascular disease, or fitted with a pacemaker, are prohibited from using them this device.
2. Patients with haemorrhage or possible haemorrhage are prohibited from using them.
3. Patients with infectious diseases are prohibited from using them this device.
4. Children and paralytics are prohibited from using them this device.
5. Those who are not sensitive to heat are prohibited from using them this device.
6. Women who are pregnant should exercise caution whilst using this device



Cleaning and Maintenance

1. These products are not waterproof, so please do not wash them. If there are stains, use soft cloth to clean it with some diluted detergent, then dry it with absorbent dry cloth.
2. Store your product in a dry and ventilated environment when you don't use it for a long time.
3. After use, please let your product cool down before folding.
4. During storage, do not put heavy things on your product.
5. Check your product frequently for damage. If damage is found, please send it to your supplier for maintenance.

Warranty Clauses

1. When used in accordance with the instructions, warranty within one year from the date of purchase.
2. Damage caused by failure to follow the instructions, human factors, force majeure, etc., is not covered by warranty.
3. Maintenance cost will be charged for product out of warranty.



Give your family a warm and healthy life.

Disclaimer Notices

In the process of the preparation of this manual, we have sought to make the content correct and complete. However, there is no guarantee that there are no errors or omissions in this manual.

We reserve the right to change specifications of the products mentioned in this manual at any time without prior declaration. No part of this manual shall be reproduced, transmitted, transcribed or stored in a searchable system or translated into any other languages without the written authorization of our company. The company is not responsible for the loss caused by non-standard or wrong operation.

Individual results may vary.

No individual result should be considered typical or guaranteed.

This product has not been evaluated by the FDA .

Leg-O-Sage System is not intended to diagnose, treat, cure, or prevent any disease, illness, or medical condition.

THIS IS NOT A MEDICAL DEVICE. This unit simply massages legs, thighs, knees, feet and arms. No medical claims are implied or warranted by use of this product. It is strongly advised that you consult your physician before use if you are currently undertaking medical treatment.

Customer Service:
Phone: 516-723-9947
Email: sales.igiany@gmail.com
www.igia.com